



**ANOGEN** - A Division of YES Biotech Laboratories Ltd.

2355 Derry Road East, Unit 23, Mississauga, ON L5S 1V6 • Tel: (905) 677-9221 • Fax: (905) 677-0023

**Mouse anti- human FPRL1 monoclonal antibody [3D7] Datasheet**

<b>Catalogue No</b>	MO-P40007A
<b>Product name</b>	Mouse anti-human formyl peptide-receptor-like-1 (FPRL1) monoclonal antibody [3D7]
<b>Clone No</b>	3D7
<b>Product type</b>	Primary antibody
<b>Description</b>	Mouse monoclonal antibody (MAb) against the residues of the 2nd extracellular loop of human FPRL1, Protein G affinity purified
<b>Target protein</b>	Human formyl peptide-receptor-like-1 (FPRL1)
<b>Immunogen</b>	The amino acid sequence of the 2nd extracellular loop of human FPRL1, conjugated with KLH.
<b>Specificity</b>	This antibody recognizes the amino acid sequence of the 2nd extracellular loop of human FPRL1 and intact human FPRL1.
<b>Reactivity</b>	Human
<b>Clonality</b>	Monoclonal
<b>Source</b>	Mouse
<b>Isotype</b>	IgG1, k
<b>Formulation</b>	Lyophilized from a solution in 0.01M PBS pH7.2
<b>Reconstitution</b>	Double distilled water is recommended to adjust the final concentration to 1.00 mg/ml. Avoid repeated freeze and thaw cycle.
<b>Applications</b>	Western blotting.
<b>Research</b>	Cellular signaling regulation and Alzheimers Disease.
<b>Storage</b>	Store at -20° C
<b>References</b>	If research is published using this product, please inform Anogen in order to cite the reference on this datasheet. Anogen will provide any product of choice as gratitude.

**This product is for LABORATORY RESEARCH USE and further manufacture only, and can not be administrated to human and animals for use in diagnostic and therapeutic procedures.**


**Manufactured by ANOGEN - A Division of YES Biotech Laboratories Limited**

**S7.5(01)**



**ANOGEN** - A Division of YES Biotech Laboratories Ltd.  
2355 Derry Road East, Unit 23, Mississauga, ON L5S 1V6 • Tel: (905) 677-9221 • Fax: (905) 677-0023

Western blot using anti-FPRL1 antibody clone 3D7:

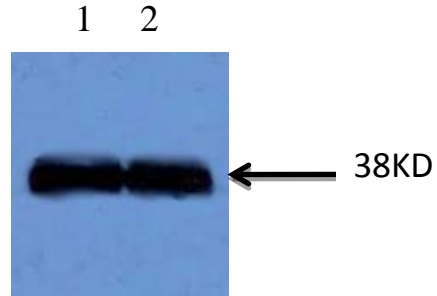

1&2: FPRL-1 /293 cells

Western blot analysis of extracts from FPRL1/293 cells using anti-FPRL1 antibody (3D7) at 1:1000 (1µg/ml) dilution.

**This product is for LABORATORY RESEARCH USE and further manufacture only, and can not be administrated to human and animals for use in diagnostic and therapeutic procedures.**

**Manufactured by ANOGEN - A Division of YES Biotech Laboratories Limited**

S7.5(01)