



ANOGEN - A Division of YES Biotech Laboratories Ltd.
2355 Derry Road East, Unit 23, Mississauga, ON L5S 1V6 • Tel: (905) 677-9221 • Fax: (905) 677-0023

Rabbit anti-Human Fibrinogen Polyclonal Antibody Datasheet

Catalogue No	AI80002
Product name	Purified rabbit anti-human Fibrinogen polyclonal antibody
Product type	Primary antibodies
Description	Rabbit anti-human Fibrinogen polyclonal antibody, purified with a Protein G affinity column
Target protein	Human Fibrinogen
Immunogen	Purified human Fibrinogen
Specificity	This antibody recognizes human Fibrinogen
Reactivity	Human, others not tested
Clonality	Polyclonal
Source	Rabbit
Formulation	Lyophilized from a solution of 0.01M PBS, pH 7.4
Applications	Suggested application and working condition: Western blot: 1:50000 dilution (antigen: 5 µg/Lane) ; ELISA 200,000~1:400,000 dilution (coating antigen: 800 ng/well) refer to the enclosed data and Experimental Fig 1. For immune-depletion, it can be used for Fibrinogen removal from test sample using affinity method.
Storage	Store in aliquots at -20°C. Avoid repeated freeze/thaw cycles
Research	Biochemistry.
References	If research is published using this product, please inform Anogen in order to cite the reference on this datasheet. Anogen will provide any product of choice as gratitude.

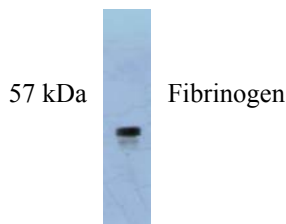


Fig.1: Western blot analysis of human fibrinogen (1.5 µg) using anti-human fibrinogen polyclonal antibody

This product is for LABORATORY RESEARCH USE and further manufacture only, and can not be administrated to human and animals for use in diagnostic and therapeutic procedures.

Manufactured by ANOGEN - A Division of YES Biotech Laboratories Limited